

# WimAdipose

## Specifications for a successful analysis

In order to achieve the best results when using WimAdipose, these are some guidelines that you should follow while conducting your assay and acquiring images from it. If your images deviate from these guidelines, the accuracy of your WimAdipose analysis measurements might be negatively affected. Although we work our best to make our solution as adaptable as possible, we do not take responsibility of their accuracy if that is the case.

Please take care of your images, we will take care of the rest.

### 1. Before acquiring your images

**Valid microscopy modalities:** WimAdipose only works with images from histological cuts. If you want to analyze images from other modalities, you can ask for a custom solution [here](#).

**Avoid artifacts:** any artifacts visible in the image (air bubbles, debris, stains...) may make it harder for the analysis to recognize your adipocytes and affect negatively your results. Please take care to avoid them when preparing your assay.

### 2. When acquiring your images


**No visible containers:** please avoid images containing parts or whole assay containers (plates, wells...). These structures can lead to false positives.

**Images in focus:** images should be properly focused in order to allow adipocytes to be viewed with the best contrast and distinguish them from the background clearly.

**Homogeneous illumination:** achieving homogeneous lightning throughout the entire image will ensure better results. Please avoid images in which some parts of the image are darker or brighter than others.

### 3. When saving your images

**Valid formats:** WimAdipose accepts images saved in jpg, jpeg, jp2, png, gif, tiff, tif, bmp, zvi, ndpi, lsm, lif and stk formats. The image uploader will not let you upload images from other formats. If you want to analyze images from a special image format, you can ask for a custom solution [here](#).

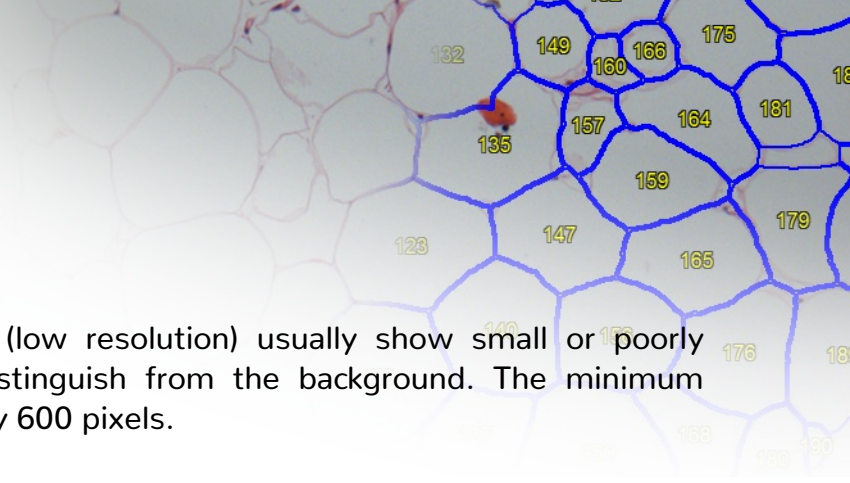

**Good image resolution:** small images (low resolution) usually show small or poorly defined objects that are difficult to distinguish from the background. The minimum resolution that we accept is 800 pixels by 600 pixels.

**No additional information visible on the image:** some microscope software add information to the images like a scale bar or a time-stamp, which can make XXX around them more difficult to detect. Please send your images without any extra information on them.